

293T cyno GITR Cell Line

Cat. No: KC-0256

Version 19022601

I.	Cell Line Information.....	1
II.	Background.....	1
III.	Cell Line Generation.....	2
IV.	Characterization using FACS.....	2
V.	Application.....	2
VI.	Cell Resuscitation.....	2
VII.	Cell Freezing.....	2
VIII.	References.....	3

I. Cell Line Information

Catalog number	KC-0256
Cell line name:	293T cyno GITR Cell Line
Gene ID/Accession #:	XM_005545123.2
Host cell line	293T
Cell type:	Human embryonic kidney
Description:	HEK293T cell line stable expressing exogenous Cyno GITR gene
Quantity:	One vial of frozen cells (5X10 ⁶ per vial)
Stability:	Stable in culture over a minimum of 10 passages
Application:	Drug screening and biological assays
Freeze medium:	70% DMEM + 20% FBS + 10% DMSO
Propagation medium:	DMEM + 10% FBS + 0.5ug/ml Puromycin
Selection marker:	Puromycin
Morphology:	Fibroblastoid cells growing as monolayer
Subculture:	Split saturated culture 1:4~1:5 every 2~3 days; seed out at about 1-3 x 10 ⁵ cells/ml
Incubation:	37 °C with 5% CO ₂
Doubling time:	Approximately 30 hours
Mycoplasma status:	Negative
Biosafety level:	1
Storage:	Liquid nitrogen immediately upon receiving

II. Background

Glucocorticoid-induced TNFR-related protein (GITR), also named Tumor necrosis factor receptor superfamily member 18 (TNFRSF18) or activation-inducible TNFR family receptor (AITR), is a cell membrane protein receptor acting as a co-stimulatory immune checkpoint molecule on T cell activation, and has proved to play a critical role in inhibiting the suppressive activity of T-regulatory cells and regulating the TCR mediated T cell death.

III. Cell Line Generation

293T Cyno G1TR cell line was generated using lentiviral vector expressing human cyno G1TR sequence.

IV. Characterization using FACS

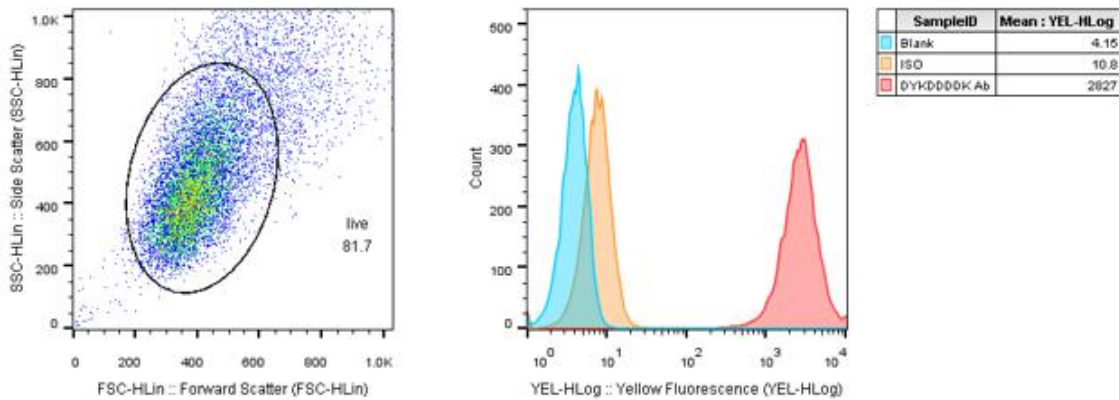


Figure: Characterization of Cyno G1TR overexpressing in 293T stable clones using ICFC FACS.

V. Application

Hybridoma or Binders of ligand screening with FACS.

VI. Cell Resuscitation

1. Prewarm culture medium (DMEM supplemented with 10% FBS and 0.5ug/ml puromycin) in a 37°C water bath.
2. Thaw the frozen vial in a 37°C water bath for 1-2 minutes.
3. Transfer the vial into biosafety cabinet and wipe the surface with 70% ethanol.
4. Unscrew the top of the vial and transfer the cell suspension gently into a sterile centrifuge tube containing 9.0 mL complete culture medium.
5. Spin at ~ 125 x g for 5~7 minutes at room temperature and discard the supernatant without disturbing the pellet.
6. Resuspend cell pellet with the appropriate volume of complete medium and transfer the cell suspension into a T25 culture flask.
7. Incubate the flask at 37°C, 5% CO₂ incubator.
8. Split saturated culture 1:4 ~ 1:5 every 2~3 days; seed out at about 1-3 x 10⁵ cells/ml.

VII. Cell Freezing

1. Prepare the freezing medium (70% DMEM + 20% FBS + 10% DMSO) fresh immediately before use.
2. Keep the freezing medium on ice and label cryovials.
3. Harvest cells to a sterile, conical centrifuge tube during the logarithmic growth period and count the cells.

4. Centrifuge the cells at 250 x g for 5 minutes at room temperature and carefully aspirate off the medium.
5. Resuspend the cells at a density of at least 3×10^6 cells/ml in chilled freezing medium.
6. Aliquot 1 ml of the cell suspension into each cryovial.
7. Freeze cells in the CoolCell freezing container overnight in a -80°C freezer.
8. Transfer vials to liquid nitrogen for long-term storage.

VIII. References

1. Nocentini, G., Ronchetti, S., Cuzzocrea, S. & Riccardi, C. GITR/GITRL: More than an effector T cell co-stimulatory system. *Eur. J. Immunol.* 37, 1165–1169 (2007).
2. Esparza, E. M. & Arch, R. H. TRAF4 functions as an intermediate of GITR-induced NF- κ B activation. *CMLS, Cell. Mol. Life Sci.* 61, 3087–3092 (2004).
3. Wyzgol, A. et al. Trimer Stabilization, Oligomerization, and Antibody-Mediated Cell Surface Immobilization Improve the Activity of Soluble Trimers of CD27L, CD40L, 41BBL, and Glucocorticoid-Induced TNF Receptor Ligand. *The Journal of Immunology* 183, 1851–1861 (2009).

Kyinno Biotechnology Co., Ltd.

Yizhuang Biomedical Park, No. 88, Beijing, China
Tel: +86-10-58222702
E-mail: bd@kyinno.com
Web: www.kyinno.com

For Research Use Only